
waterpik[®] SinuSense™


QUICK START GUIDE

WATER PULSATOR


ENGLISH

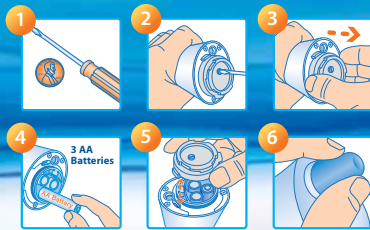1. To install batteries: locate the gray battery cover on the bottom of the white handle.
2. Use a standard screwdriver or coin to loosen (but not remove) the center screw on the gray battery cover.
3. It will take approximately five complete turns. Note: Screw will remain in battery cover. Pull the gray battery cover off the white handle to reveal the battery compartment.
4. Install three AA batteries following the diagram inside the battery compartment.
5. Reinstall the gray battery cover by aligning the tab on the gray battery cover with the slot on the white handle, and tighten screw.
6. Install SoftSeal® nozzle onto nozzle tip. See instruction manual for details.

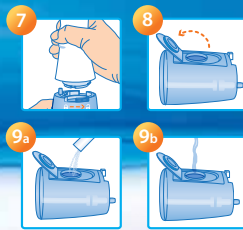7. Remove the blue reservoir (if attached to the white handle) by sliding the lock on the top of the blue reservoir to the unlocked position and pull the white handle away from the blue reservoir.
8. On a flat surface, place the blue reservoir on its side with the door facing up, then flip open the door.
9. (9a) Tear open a Soothing Saline Easy-Pour Pack and empty contents into the blue reservoir. (9b) Fill the reservoir full with 8 ounces of distilled water or previously boiled water then press the reservoir door closed.
Important: Be sure to turn reservoir upright after filling or seepage from the vent may occur.

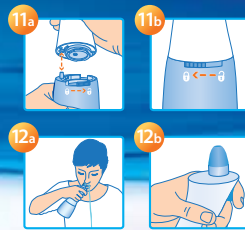10. Heat blue reservoir in microwave at 5-second intervals to desired temperature. **Do not microwave longer than 30 seconds.**
11. (11a) Align the tube on the blue reservoir with the hole in the white handle, (11b) press together and slide the lock on top of the blue reservoir to the locked position. Gently shake to mix the solution.
12. (12a) While at the sink, lean forward and fit the nozzle of the water pulsator inside the right nostril. (12b) Control the water flow by squeezing gently on the FlowControl® trigger for 10 seconds (or until 4 ounces of solution is used). To stop the flow, release the trigger. Repeat as needed on both sides.
13. After use, gently shake the blue reservoir to discard the remaining solution. Leave reservoir door open to air dry.

Helpful Tips

- Remember to breathe through your mouth
- Gently blow your nose after rinsing both sides
- Use daily for best results

Go to

www.waterpiksinusense.com
to learn more about how to use Waterpik® SinuSense™ products or refer to your device instruction manual for more information.

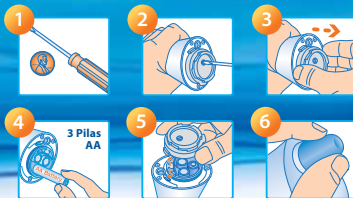
**Register your product online now!
Receive up to \$42 in savings on the
purchase of additional Waterpik®
products.**

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© 2012 Water Pik, Inc.