

waterpik®
waterflosser™

Congratulations!

You are now on your way
to better oral health!

✓ Register and Save

Register your product at
*waterpik.com/registration to receive over **\$50**
in savings (includes **50%** off your next order of
assorted tips + tip storage case)

*For Canadian customers, register at
waterpik.ca/registration.

✓ Write a Review

Creating innovative Waterpik® products
that inspire better oral health is our
passion. We invite you to *leave a review by
sharing your thoughts on your new Water
Flosser on waterpik.com or the retailer's
website.

*US customers only

✓ We're Here to Help

We hope you are loving your new Water Flosser.
If for any reason you are not satisfied or have
any product questions, please contact our fast
and friendly Customer Service Team based in
Fort Collins, Colorado. Your new Water Flosser
is backed by a 1-year manufacturer's warranty.

Call: 1-800-525-2774 (USA)
1-888-226-3042 (CAN)
Monday-Friday
7 am - 5 pm MST

Chat or E-mail
service@waterpik.com



**Quick Start
Guide Inside!**



Register and Save!

Save over \$50!

Register at ***Waterpik.com/Registration** to
receive exclusive discount code.

*For Canadian customers, register at
waterpik.ca/registration

Offer includes
50% off
your next order
of assorted
tips + tip
storage case.



We'd Love to Hear from You!



Connect with Us
@[Waterpik](https://www.instagram.com/waterpik)



Call: 1-800-525-2774 (USA) 1-888-226-3042 (CAN)
www.waterpik.com www.waterpik.ca

Waterpik® 14 Day Challenge App




Make water flossing a habit! Try our new mobile App
for the best water flossing experience.




Quick Start Guide

GETTING STARTED

- 1 Remove reservoir.

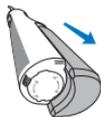

- 2 To remove battery door, twist, and pull.

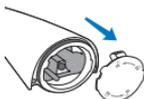

- 3 Insert 3 AA batteries.



- 4 To replace battery door, line-up arrows, push and twist.



- 5 Insert tip firmly until it clicks into place.



- 6 If the reservoir was removed, insert reservoir back on the unit. Fill reservoir with warm water.



WATER FLOSS

- 7 With water switch in OFF position, place tip in mouth. Slide power switch to low (I) and increase to your liking.




- 8 Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.



Keep unit upright during use to ensure continuous water flow.

TECHNIQUE

- 9 Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

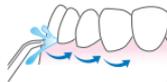

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

NEED HELP?

Visit Waterpik.com/WelcomWF
To watch a How To Video and find additional product support.

**Call Us! 1-800-525-2774 (USA)
1-888-226-3042 (CAN)**