

## QUICK START GUIDE

Review full instruction book prior to first use.

## 1 FIRST USE

1. Remove reservoir.



2. To remove battery door, twist, and pull.



3. Insert 3AA batteries. Included.



4. To replace battery door, line-up arrows, push and twist.



5. Insert tip firmly until it clicks into place.



6. Fill reservoir with warm water.



## 2 START WATER FLOSSING

**TO AVOID MESS: PLACE THE TIP IN YOUR MOUTH BEFORE TURNING UNIT ON**

Slide power switch to low(I) or high(II) position. Start with lowest setting and increase to your liking.



## 3 LEAN OVER SINK

Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.



## 4 TECHNIQUE

Aim water just above the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.



## 5 HOW TO VIDEO

Visit [waterpik.com/vid02](http://waterpik.com/vid02)

Register your product and receive over \$50 in savings.  
[waterpik.com/registration](http://waterpik.com/registration)

# GUÍA RÁPIDA PARA COMENZAR

Revise el manual de instrucciones por completo antes del primer uso.

## 1 PRIMER USO

1. Quite el depósito de agua.



2. Para quitar la puerta de la batería, gírela y jale.



3. Introduzca 3 pilas AA (incluidas).



4. Para volver a colocar la tapa de las pilas, alinee las flechas, empuje y gire.



5. Inserte el cabezal con firmeza hasta que encaje en su lugar.



6. Llene el depósito con agua tibia.



## 2 INICIE EL IRRIGADOR BUCAL


**EVITE DERRAMES:  
COLOQUE EL CABEZAL EN  
LA BOCA ANTES DE  
ENCENDER LA UNIDAD**

Deslice el interruptor de encendido hacia la posición baja(I) o alta(II). Inicie con el ajuste más bajo e increméntelo según su preferencia.



## 3 INCLÍNESE SOBRE EL LAVABO

Cierre los labios lo suficiente para evitar salpicar mientras el agua continúa fluyendo de la boca hacia el lavabo.

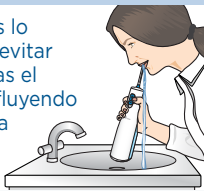

## 4 TÉCNICA

Apunte el agua justo por encima de las encías en un ángulo de 90 grados. Siga la línea de las encías y haga pausas breves entre los dientes.



## 5 CÓMO VER VIDEOS

Visite [waterpik.com/vid02](http://waterpik.com/vid02)

Registre su producto y reciba más de USD 50 en descuentos.  
[waterpik.com/registration](http://waterpik.com/registration)