A dynamic splash of clear blue water with numerous bubbles, moving from left to right across the middle of the page. A solid magenta circle is partially visible on the right edge of the page.

# waterpik®

## SIDEKICK™

Instructions For Use  
Instructions d'utilisation  
Instrucciones de uso

**We are here for you | Nous sommes à votre disposition | Estamos a su disposición**  
1-800-525-2774 (USA) 1-888-226-3042 (Canada) [waterpik.com/HowTo](http://waterpik.com/HowTo)

## IMPORTANT SAFEGUARDS

---

### IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

#### READ ALL INSTRUCTIONS BEFORE USING.

#### DANGER:

##### To reduce the risk of electrocution:

- Do not handle plug with wet hands.
- Do not place in or drop into water or other liquid.
- Do not use while bathing.
- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not reach for a product that has fallen into water. Unplug immediately.
- Check the power supply and the power supply cord for damage before the first use and during the life of the product.

### WARNING:

#### To reduce the risk of burns, electrocution, fire, or injury to persons:


- If product is opened/disassembled for ANY reason, warranty is **VOID**.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Water Pik, Inc. or your local distributor which can be found at [www.waterpik.com](http://www.waterpik.com).
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
- Use this product only as indicated in these instructions or as recommended by your dental professional.
- Only use tips and accessories recommended by Water Pik, Inc.
- Do not drop or insert any foreign object into any opening or hose.
- Keep electrical cord away from heated surfaces.

- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not use iodine or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product. Use of these products may void warranty.
- Fill reservoir with water or other dental professional recommended solution only.
- Remove any oral jewelry prior to use of this product.
- Do not use if you have an open wound on your tongue or in your mouth.
- If your physician or cardiologist has advised you to receive antibiotic premedication before dental procedures, you should consult your dentist before using this instrument or any other oral hygiene aid.
- Close supervision is necessary when this product is used by, on, or near children or invalids.
- Do not use this product for more than 3 minutes in each two-hour period.
- Unit will shut off automatically after 3 minutes of continuous usage.
- Children should be supervised to ensure that they do not play with the appliance.
- Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet.


### **SAVE THESE INSTRUCTIONS**

## PRODUCT DESCRIPTION

1. Reservoir
2. Base unit
3. Water hose release button
4. Water hose connection
5. Water hose
6. Reservoir valve
7. Handle
8. Handle notch
9. Power button
10. Pressure control
11. Power supply storage
12. Power supply
13. Power supply cord
14. Travel case

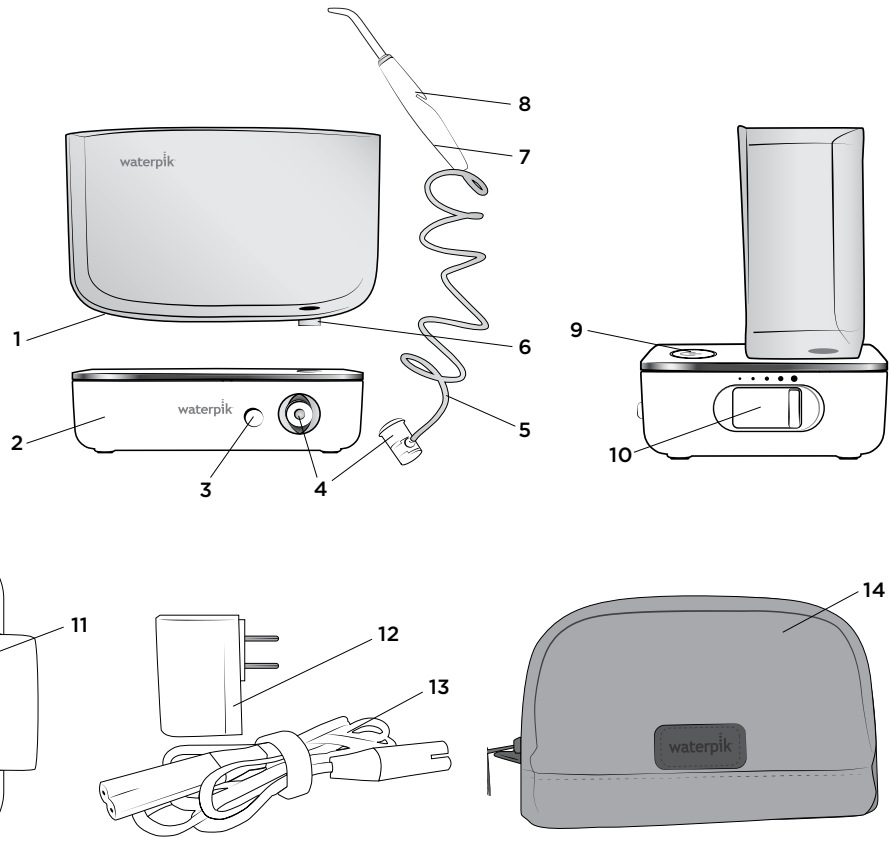Replacement accessories may be purchased online at [www.waterpik.com](http://www.waterpik.com), or by phone from Water Pik, Inc. at 1-800-525-2774 (USA) or 1-888-226-3042 (Canada). To find detailed information about tips, other accessories, and product usage please visit our website at [www.waterpik.com](http://www.waterpik.com).
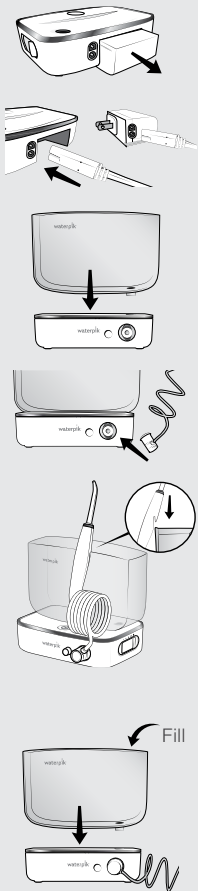
### Assembly



1. Remove power supply from power supply storage.
2. Connect power supply cord to power supply and to base unit.
3. Place reservoir on base unit and press down firmly until secured.
4. Connect water hose to base unit.
5. Using handle notch, hook handle onto front edge of reservoir.

### Powering the Unit

Plug the power supply into an electrical outlet. If the outlet is controlled by a wall switch, be sure it is turned on.

### Preparing the Reservoir

Remove reservoir from base unit and fill with warm water. Replace reservoir onto the base unit and press down firmly.
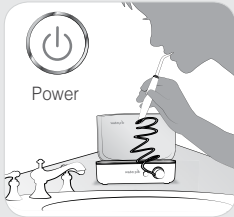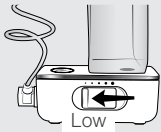
### Prime Unit (First Time Use Only)

Set pressure control switch to high, point the handle and tip into sink.  
Turn power ON until water flows from the tip. Turn power OFF.

### Adjusting the Pressure Setting



Slide the pressure control on the base of the unit to the lowest setting for first-time use.  
Gradually increase pressure over time to the setting you prefer or as instructed by your dental professional.

### Using the Water Flosser

Always place the tip in your mouth before turning the water ON.  
Lean over sink and slightly close lips enough to prevent splashing while still allowing water to flow from mouth into the sink. Turn power ON. When finished, turn power OFF.

### Recommended Technique

For best results, start with back teeth and work toward the front teeth.  
Direct the stream at a 90-degree angle to your gumline with tip in mouth.  
Glide tip along gumline and pause briefly between teeth.  
Continue until you have cleaned the inside and outside of both the upper and lower teeth.
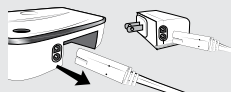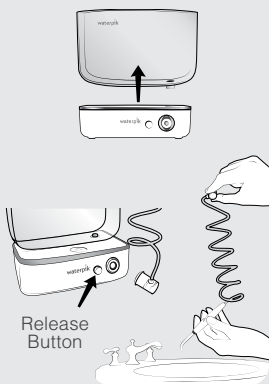
### Disassembly for Travel

1. Remove and empty reservoir.
2. Unplug unit from electrical outlet.
3. Disconnect handle hose from base unit by pressing the water hose release button.  
Point tip into sink and extend water hose to drain all remaining water.
4. Disconnect power supply cord from base unit and power supply.
5. Insert power supply into power supply storage.
6. Slide base unit into reservoir.

### Travel

Place the full unit and all components inside the provided travel case for protection when traveling.

## LIMITED THREE-YEAR WARRANTY

---

### **Service Maintenance**

Waterpik® Water Flossers have no consumer serviceable electrical items and do not require routine service maintenance. For all parts and accessories go to [www.waterpik.com](http://www.waterpik.com). You may also call our toll-free customer service number at 1-800-525-2774 (USA) or 1-888-226-3042 (Canada).

Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Still have questions? Go to [www.waterpik.com](http://www.waterpik.com).

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 3 years from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, which in its opinion is defective, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Installation is the responsibility of the consumer and is not covered by the warranty. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

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**[www.waterpik.com](http://www.waterpik.com)**

Waterpik® is a trademark of Water Pik, Inc. registered in Argentina, Australia, Benelux, Brazil, Canada, Chile, China, Czech Republic, Denmark, EU, Finland, France, Germany, Hong Kong, Hungary, India, Israel, Italy, Japan, Kazakhstan, Korea, Malaysia, Mexico, New Zealand, Norway, Pakistan, Philippines, Poland, Russian Federation, Saudi Arabia, Singapore, South Africa, Spain, Sweden, Switzerland, Taiwan, Turkey, Ukraine, United Kingdom, United States and Vietnam.

Waterpik® (stylized) is a trademark of Water Pik, Inc. registered in Australia, Brazil, Canada, China, EU, Hong Kong, India, Israel, Japan, Korea, Malaysia, Mexico, Philippines, Russian Federation, Saudi Arabia, South Africa, Switzerland, Taiwan, Thailand, Trinidad & Tobago, Ukraine, United States, and Vietnam.

SIDEKICK™ is a trademark of Water Pik, Inc.

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SIDEKICK™ est une marque de commerce de Water Pik, Inc.
SIDEKICK™ es una marca de comercio de Water Pik, Inc.