


waterpik® waterflosser®

QUICK START GUIDE WP-100 Series


Register your product and receive over \$50 in savings.
waterpik.com/registration

Visit waterpik.com/welcomeWF
to watch How To Video

1 Unit Prep

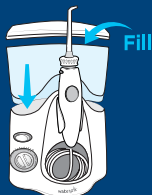

A. Remove the reservoir from base and fill with lukewarm water. Replace the reservoir firmly on the base.

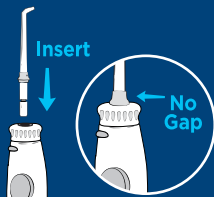


B. Press tip firmly into handle. There should be no gap between handle and colored tip ring.



C. Turn pressure to high, point tip into sink, turn unit **ON** until water flows. Turn unit **OFF**.


2 Water Floss

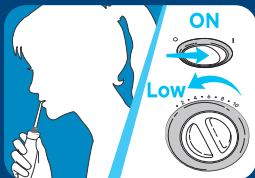

A. Start on low pressure and increase to your liking. With unit **OFF**, place tip in mouth. Turn power **ON**.



B. Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.

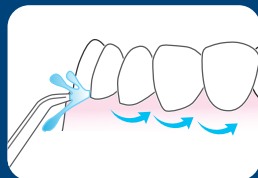

C. Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

www.waterpik.com

This device is only compatible with voltages of 120 VAC, 60Hz

©2018 Water Pik, Inc.

waterpik® waterflosser®

GUÍA RÁPIDA PARA COMENZAR WP-100 Serie

Registre su producto y reciba
más de USD 50 en descuentos.
waterpik.com/registration

Visite waterpik.com/welcomeWF
para ver el video práctico

1 Preparación de la unidad



A. Retire el depósito de la base y llene con agua tibia. Vuelva a colocar el depósito firmemente sobre la base.




B. Presione con firmeza la punta dentro del mango. No debería haber espacio entre el mango y el aro de punta de color.



C. Gire la presión hacia Alto, dirija la punta hacia el lavabo, encienda la unidad hasta que el agua fluya. Apague la unidad.


2 Cómo utilizar el Irrigador oral

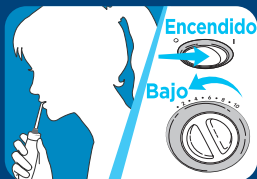

A. Comience con presión baja e increméntela según su agrado. Con la unidad apagada, coloque la punta en la boca. Encienda la unidad.



B. Cierre los labios lo suficiente para evitar salpicar mientras el agua continúa fluyendo de la boca hacia el lavabo.

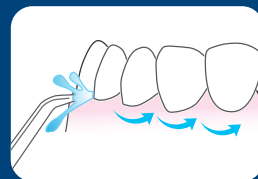

C. Dirija el agua hacia las encías en un ángulo de 90 grados. Siga las encías y haga pausas breves entre los dientes.

EVITE DERRAMES: MANTENGA LA PUNTA EN LA BOCA MIENTRAS LA UNIDAD ESTÁ ENCENDIDA.