

# waterpik® waterflosser

Model WP-260

## QUICK START GUIDE

### 1 Prepare the Reservoir:



Remove the reservoir from base and fill with lukewarm water.

Replace



Replace the reservoir onto the base, **PRESSING DOWN FIRMLY** to ensure full seating.

### 2 Select and Insert Tip:




Press tip firmly into handle.



When tip is properly inserted there should be no gap between handle and colored tip ring.


### 3 Prime Unit (First Time Use Only):

**Prior to Very First Use of Unit, Prime Unit:** Set pressure control dial on base to **high**, point the handle and tip into sink. Turn unit **ON** until water flows from the tip. Turn unit **OFF**.

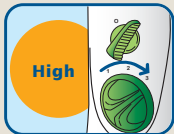

### 4 Adjust Pressure Control:

For first time users, adjust pressure control to **lowest position**. Gradually increase pressure to your liking.

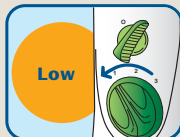

*continued on back*


This device is compatible with global voltages of 100 - 240VAC, 60/50Hz.

5


## Using your Waterpik® Water Flosser:

**Always place tip in mouth before turning unit on.**

Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. Turn unit **ON**.

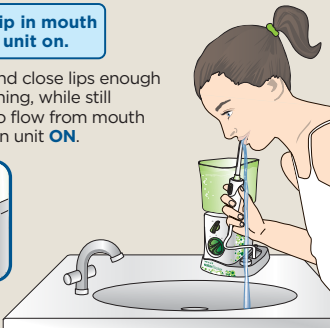

Aim the water stream just above the gumline at a 90 degree angle.

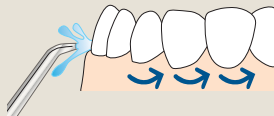

Starting with the back teeth, follow the gumline and pause briefly between teeth.

Continue until all areas around and between teeth have been cleaned. Turn unit **OFF**.

6

## Tip Removal:

Slide the tip eject button on the handle. Remove the tip.

**Caution:** Do not slide tip eject button while unit is running.



7

## REGISTER YOUR PRODUCT ONLINE NOW!

Receive up to  
\$44

in savings on the purchase of additional Waterpik® products.  
[www.waterpik.com/registration](http://www.waterpik.com/registration)