

waterpik® waterflosser

Model
WP-300

QUICK START GUIDE

1 Prepare the Reservoir:



Remove the reservoir from base and fill with lukewarm water.

Replace



Replace the reservoir onto the base, **PRESSING DOWN FIRMLY** to ensure full seating.

2 Select and Insert Tip:



Press tip firmly into handle.



When tip is properly inserted there should be no gap between handle and colored tip ring.


3 Prime Unit (First Time Use Only):

Prior to Very First Use of Unit, Prime Unit: Set pressure control dial on base to **high**, point the handle and tip into sink. Turn unit **ON** until water flows from the tip. Turn unit **OFF**.



4 Adjust Pressure Control:

For first time users, adjust pressure control to **lowest position**. Gradually increase pressure to your liking.

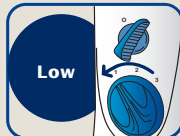

continued on back


This device is compatible with global voltages of 100 - 240VAC, 60/50Hz.

5


Using your Waterpik® Water Flosser:

Always place tip in mouth before turning unit on.

Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. Turn unit **ON**.

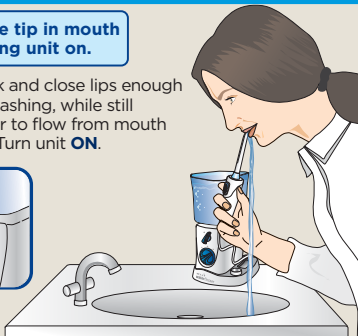

Aim the water stream just above the gumline at a 90 degree angle.

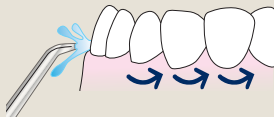

Starting with the back teeth, follow the gumline and pause briefly between teeth.

Continue until all areas around and between teeth have been cleaned. Turn unit **OFF**.

6

Tip Removal:

Slide the tip eject button on the handle. Remove the tip.

Caution: Do not slide tip eject button while unit is running.



7

REGISTER YOUR PRODUCT ONLINE NOW!

Receive up to
\$44

in savings on the purchase of additional Waterpik® products.

www.waterpik.com/registration