

waterpik® waterflosser®


QUICK START GUIDE

WP-560/562/563/564/565/566/569


Register your product and receive over \$50 in savings.
waterpik.com/registration

Visit waterpik.com/vid560
to watch How To Video

1 Unit Prep

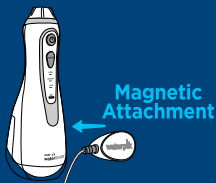


A. Charge unit for 4 hours by placing the magnetic charger on the front of the unit.

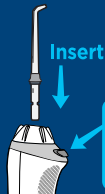


Caution:
Do not press tip eject button while unit is running.

B. Insert tip until it clicks into place. To remove tip, push the tip eject button.

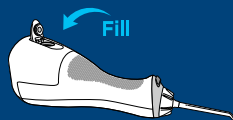

C. Lift the flip top on reservoir and fill with warm water. Refill as needed.


2 Water Floss

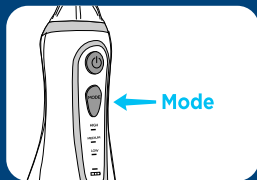

A. Press mode button to adjust pressure. Start on low pressure and increase to your liking.



B. Turn power ON. Close lips enough to prevent splashing. Keep unit upright during use to ensure continuous water flow.

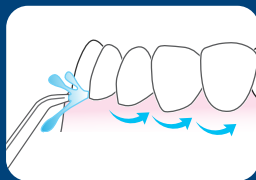

C. Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

waterpik waterflosser®

GUIDE DE DÉMARRAGE RAPIDE

WP-560/562/563/564/565/566/569

Enregistrez votre produit
et bénéficiez de plus de
50 \$ d'économies.
waterpik.com/registration


Veuillez visiter le site Web
waterpik.com/vid560 pour
visionner la vidéo d'instructions

1 Préparation de l'appareil



Fixation
magnétique


A. Chargez l'appareil pendant 4 heures en plaçant le chargeur magnétique sur l'avant de l'appareil.

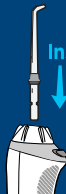

Insertion

Mise en garde :
N'appuyez pas sur le
bouton d'éjection de
l'embout lorsque
l'appareil est en marche.

B. Insérez l'embout jusqu'à ce qu'il s'emboîte. Pour retirer l'embout, appuyez sur le bouton d'éjection de l'embout.

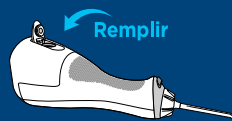

Remplir

C. Soulevez le couvercle rabattable du réservoir et remplissez le réservoir d'eau tiède. Remplissez à nouveau le réservoir lorsque cela est nécessaire.

2 Utilisation de l'hydropulseur




Mode

A. Appuyez sur le bouton de sélection du mode pour régler la pression. Démarrez sur le réglage le plus bas et augmentez la pression selon vos préférences.



B. ALLUMEZ l'appareil. Gardez vos lèvres suffisamment fermées pour éviter les éclaboussures. Laissez l'appareil en position verticale pendant l'utilisation pour assurer un débit d'eau continu.

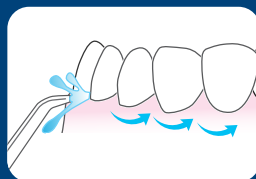

C. Dirigez le jet d'eau juste au-dessus de la ligne gingivale à un angle de 90 degrés. Suivez la ligne gingivale et faites de courtes pauses entre les dents.

POUR ÉVITER LES ÉCLABOUSSURES : GARDEZ L'EMBOU DANS LA BOUCHE LORSQUE L'APPAREIL EST EN MARCHE.