

waterpik® waterflosser®


English

QUICK START GUIDE

Review full instruction book prior to use


1 FIRST USE

Charge unit for 4 hours by placing the magnetic charger on the front of the unit.

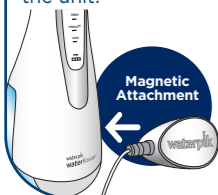

Insert tip firmly until it clicks into place.

Fill reservoir with warm water.

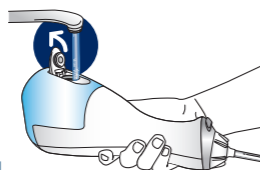

2 START WATER FLOSSING

Press mode button to adjust pressure. Start with the lowest setting and increase to your liking.

**TO AVOID MESS:
PLACE THE TIP
IN YOUR MOUTH
BEFORE TURNING
UNIT ON**



3 LEAN OVER SINK

Close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink.



4 TECHNIQUE

Aim water just above the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.



5 HOW TO VIDEO



Scan here or visit waterpik.com/OH560

Register your product and receive over \$30 in savings. waterpik.com/registration

Model WP-560/562/563/564/565/566

This device is compatible with voltages of 100-240VAC, 50/60 Hz

©2016 Water Pik, Inc.

1.800.525.2774 (USA) 1.888.226.3042 (Canada)

waterpik® waterflosser®

Français

GUIDE DE DÉMARRAGE RAPIDE

Lisez l'intégralité du manuel d'instructions avant utilisation


1 PREMIÈRE UTILISATION

Chargez l'appareil pendant 4 heures en plaçant le chargeur magnétique sur l'avant de l'appareil.



Insérez fermement l'embout jusqu'à ce qu'il s'emboîte.

Remplissez le réservoir d'eau tiède.

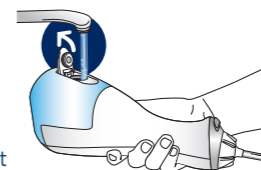

2 DÉMARRAGE DE L'HYDROPROPULSEUR

Appuyez sur le bouton de sélection du mode pour régler la pression. Démarrez sur le réglage le plus bas et augmentez la pression selon vos préférences.

**POUR ÉVITER LES ÉCLABOUSSURES :
PLACEZ TOUJOURS
L'EMBOUT DANS LA
BOUCHE AVANT DE
METTRE L'APPAREIL
EN MARCHÉ.**



3 PENCHEZ-VOUS PAR-DESSUS L'ÉVIER

Gardez vos lèvres suffisamment fermées pour empêcher les éclaboussures, tout en laissant l'eau s'écouler de votre bouche dans le lavabo.



4 TECHNIQUE

Dirigez le jet d'eau juste au-dessus de la ligne gingivale à un angle de 90 degrés. Suivez la ligne gingivale et faites de courtes pauses entre les dents.



5 VIDÉO D'INSTRUCTIONS



Scannez ce code ou visitez le site Internet waterpik.com/OH560

Enregistrez votre produit et recevez plus de 30 \$ d'économies. waterpik.com/registration

Modèle WP-560/562/563/564/565/566

Cet appareil est compatible avec les valeurs de tension de 100-240VAC, 50/60 Hz

20019833-F AF

1.800.525.2774 (É.-U.) 1.888.226.3042 (Canada)