



QUICK START GUIDE

1 Prepare the Reservoir:

Lift the hinged reservoir lid and remove the reservoir from base and fill with lukewarm water.

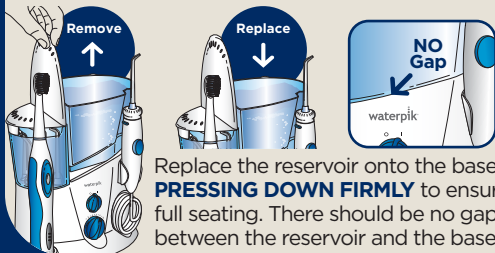


2 Select and Insert Tip:

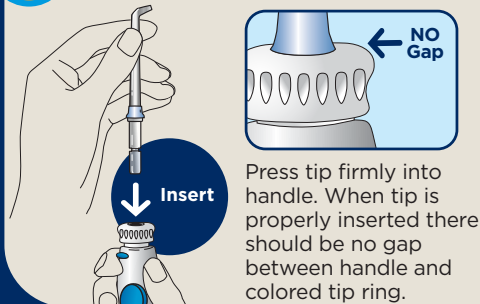

3 Prime Unit (First Time Use Only):

Prior to Very First Use of Unit, Prime Unit:


Set pressure control dial on base to **high**, point the handle and tip into sink. Turn unit **ON** until water flows from the tip, turn unit **OFF**.

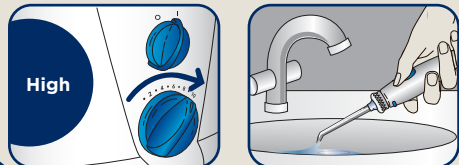

4 Adjust Pressure and Pause Control:

For first time users, adjust pressure control to **lowest position**. Gradually increase pressure to your liking.

You can temporarily stop water flow by pressing the pause button on the handle.

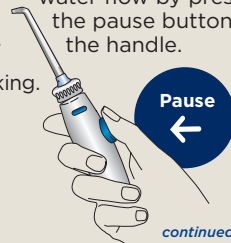


continued on back


5 Using your Waterpik® Water Flosser:

Always place tip in mouth before turning unit on.

Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. Turn unit **ON**.

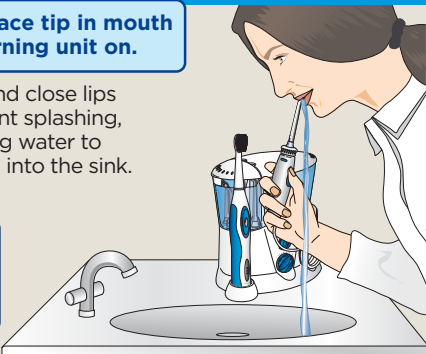

Aim the water stream just above the gumline at a 90 degree angle.

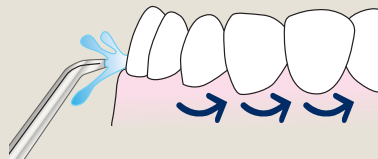



Starting with the back teeth, follow the gumline and pause briefly between teeth.

Continue until all areas around and between teeth have been cleaned.

6 Tip Removal:

Push the tip eject button on the handle. The control knob will rise, releasing the tip. Remove the tip.

Caution: Do not press tip eject button while unit is running.

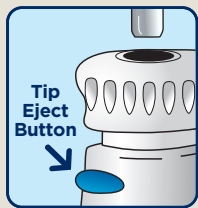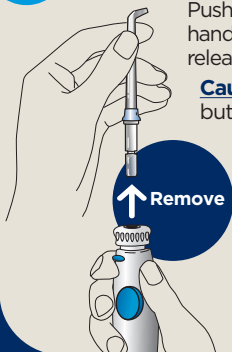

7 REGISTER YOUR PRODUCT ONLINE NOW!

Receive up to
\$44

in savings on the purchase of additional Waterpik® products.
www.waterpik.com/registration

Water Pik, Inc. 1.800.525.2774 • www.waterpik.com

1730 East Prospect Road, Fort Collins, CO 80553-0001 USA

©2012 Water Pik, Inc.

20016219-F AA