

waterpik® waterflosser®


QUICK START GUIDE

Model WP-940

Register your product and receive over \$50 in savings.
<http://www.waterpik.com/orthocare/registration>

For full instructions, review instruction manual or visit
waterpik.com/OH940 and watch the How To Video

Step 1 Unit Prep

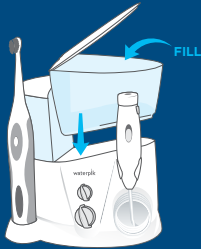

A. Fill reservoir with warm water and replace it firmly on the base.



B. Insert tip firmly until it clicks into place.



C. Turn pressure to high, point tip into sink, turn power and water switches on until water flows.

Step 2 First Use



A. Place tip in mouth and turn power on.




B. With tip in mouth, turn handle switch to ON position.



C. Start on low pressure and increase to your liking.


TO AVOID MESS: KEEP TIP IN MOUTH WHILE UNIT IS ON.

Step 3 Treatment Phase: Recommended Tip Use

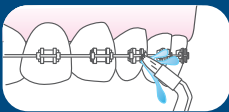

Orthodontic Tip

Lightly brush area between teeth and all around orthodontic bracket.

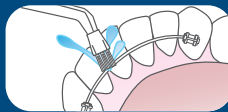

Plaque Seeker® Tip

Place the tip close to the teeth so that the bristles are gently touching the teeth.



Classic Jet Tip

Aim water just above the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.